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**User's Guide**

Advanced Sonic Pulse Electric Toothbrush

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**AQ-110**

## CONTENTS

Important Safe Guards	4
Configuration	10
Getting Started	12
How To Use	16
Sanitizing the Brush Head	18
Maintenance	20
Trouble Shooting	22
Disposal	24

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## IMPORTANT SAFE GUARDS

The Power toothbrush is safe and effective when used according to directions. Like any electrical appliance, however, basic safety precautions should be followed, especially when children are present.

### READ ALL INSTRUCTIONS BEFORE USING

## 1. DANGER - WATER

- 1 Do not place or store the charger base with UV sanitizer where it can fall or be pulled into a bathtub, sink or toilet.
- 2 Do not place or drop the charger into water or other liquid.
- 3 Do not reach for a product that has fallen into water or other liquid. Unplug immediately.
- 4 Do not use while bathing.
- 5 Unplug charger base with UV sanitizer before cleaning and assure that it is dry before plugging it into an electrical outlet.

## 2. DANGER - VOLTAGE

- 1 Do not plug the charger base with UV sanitizer into an outlet with a voltage other than specified on the bottom of the charger.
- 2 This product is designed to operate within a range of 220 volts.
- 3 Voltage converters do not guarantee compatibility.
- 4 Do not force the plug into an outlet.

## 3. DANGER - WARNINGS

- 1 This product is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities or lack of experience and knowledge, unless they have been given adequate supervision or instruction concerning use of the product by an adult responsible for their safety.

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- 2) Children should be supervised to ensure they do not play with the product.
- 3) Use this product only for its intended use as described in this manual. This product is designed to clean your teeth, gums and tongue only. Do not use any attachments or brush heads other than those recommended by the manufacturer.
- 4) The toothbrush is a personal care appliance and is not intended for use on multiple patients in a dental office or institution.
- 5) If your toothpaste includes peroxide, baking soda or other bicarbonate (common in whitening toothpastes), thoroughly clean your brush head and handle to avoid cracking the plastic.
- 6) Never operate this product if it has a damaged cord, plug or charger base with UV sanitizer. Do not use this product if it is not working properly, if it has been dropped or damaged or dropped into water. Immediately discontinue use if the product appears damaged in any way.
- 7) Keep the toothbrush, charger base and power cord away from heated surfaces.
- 8) If you experience discomfort or pain in your mouth, discontinue use of this product and contact a dentist.
- 9) Do not attempt to remove the battery.
- 10) Do not clean the toothbrush in the dishwasher.
- 11) Discontinue use of the charger base with UV sanitizer if the UV light bulb remains on when the lid is open. UV light can be harmful to human eyes and skin. Keep the charger base with UV sanitizer out of the reach of children at all times.
- 12) The UV light bulb is hot during and immediately after the sanitizing cycle. Do not touch the UV light bulb when it is hot.
- 13) Do not insert brush heads or operate the UV sanitizer without the UV sanitizer drip tray firmly in place.
- 14) Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.

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- 15) Do not leave appliance when plugged in. Unplug from outlet when not in use and before servicing.
- 16) Do not pull or carry by cord, use cord as a handle, close a door on cord, or pull cord around sharp edges or corners. Do not run appliance over cord.
- 17) Do not unplug by pulling on cord. To unplug, grasp to plug, not the cord.
- 18) Do not handle plug or appliance with wet hands.
- 19) Turn off all controls before unplugging.
- 20) Replacement of the UV lamp can only be carried out by professional staff, can not be operated by the user.
- 21) This cover is provided with an interlock to reduce the risk of excessive ultraviolet radiation.

## 4. MEDICAL-WARNINGS

- 1) If you have had oral surgery or gum surgery within the past two months, consult your dental provider before using the product.
- 2) Consult your dental provider if excessive bleeding occurs while using this product or minor bleeding continues to occur after one week of use.
- 3) Consult your physician or implant device manufacturer prior to using this product if you have questions concerning the combined use of this product with a pacemaker or other implanted medical device(s).
- 4) Consult your physician prior to using the toothbrush if you have any medical concerns.
- 5) The toothbrush has been tested and is certified as per KC(Korean certification for safe electric appliances) regulation.
- 6) Never use while sleeping or drowsy.

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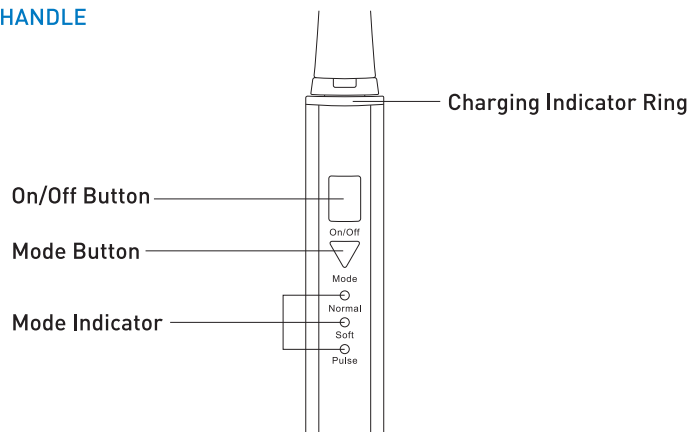
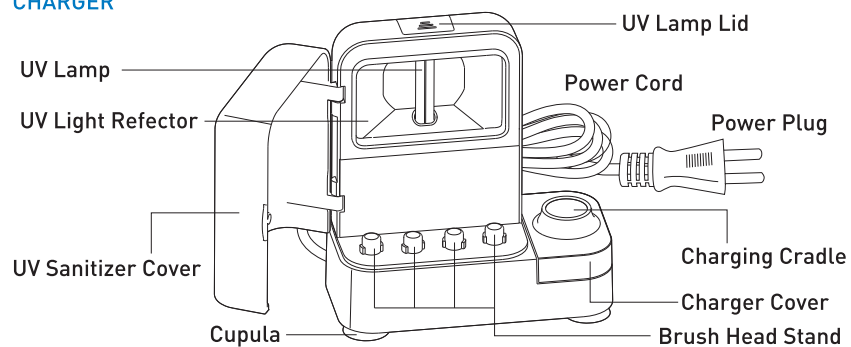
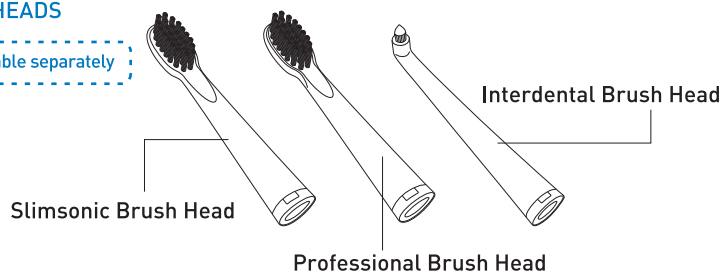
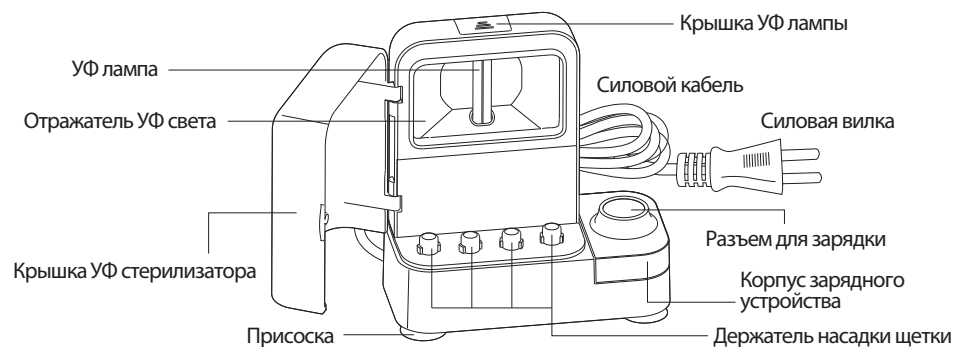
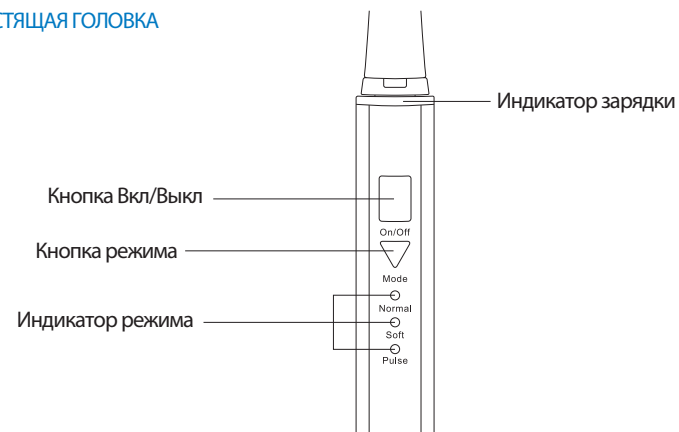
## GETTING STARTED

### Installing the charger Power Cord

Place the cord in the right, left or center cord guide, depending on where the charger will be placed.

### Place the charger on the surface

Chose a clean, level, smooth surface to place charger, and gently pressing down the cupules at bottom of the charger to suck the surface.




### Mounting on the wall

Purchase two round wood screws (diameter: 2.8~3.2mm), fasten the screw on the wall, so that they extend about 2.5mm from the wall and are positioned 74mm apart.

### Charging the toothbrush

Be sure to turn off the toothbrush before charging.

Place the handle in the charger.  
The charging indicator ring will light up red.

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
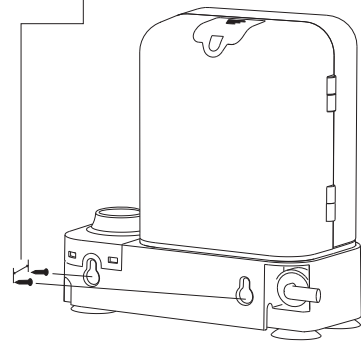
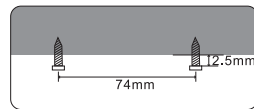
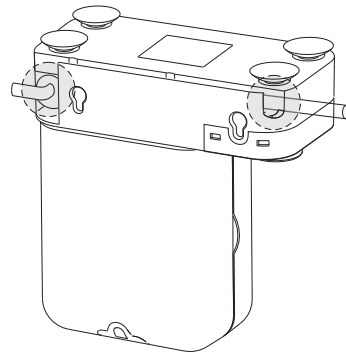
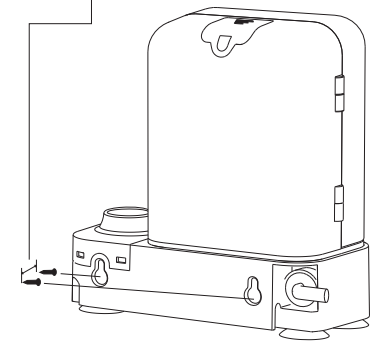
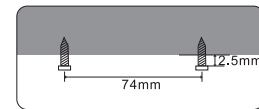
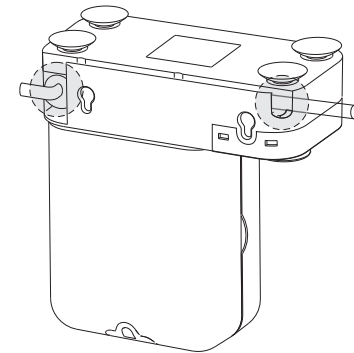
#### Notes:

- 1 It takes at least 24 hours to charge the battery fully. The indicator light will change into blue after fully charged.
- 2 The toothbrush must be charged for 24 hours, if:
  - a. this is the first time using the toothbrush;
  - b. the toothbrush has not been used for more than three months; or
  - c. if the charging indicator blinks red during use.
- 3 To keep the battery fully charged, we recommend that you keep the toothbrush in the charger when not in use.
- 4 During charging, the toothbrush will get warm, but this is normal.
- 5 The operating duration will become shorter as the battery gets older. If the toothbrush operates only for 1 short time after being fully charged, this means that the battery is nearing the end of its life cycle. Contact service center for support.

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



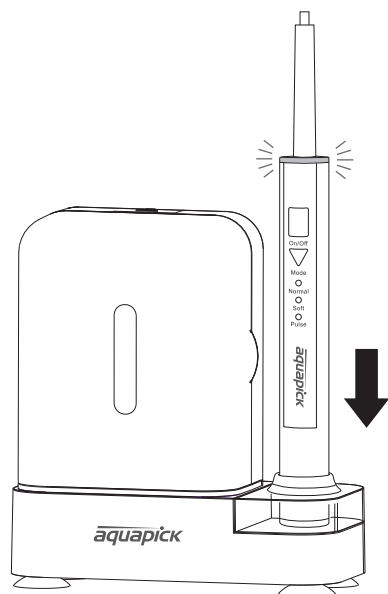
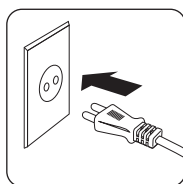
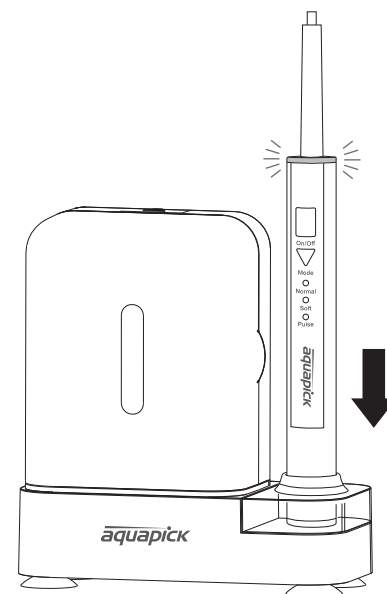
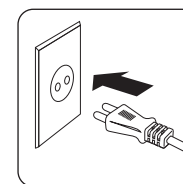
## HOW TO USE

- 1 Align the slot of brush head to the protrusion of handle and slide down, and then turn the brush head until it clicks into place.(Fig 1)
- 2 Wet the bristles and apply a small amount of toothpaste.
- 3 Place the bristles against the teeth at a slight angle towards the gum line.(Fig 2)
- 4 Press the mode button to select mode.
  - a. Normal: Standard mode for superior teeth cleaning.
  - b. Soft: Gently cleans your teeth or any other sensitive areas.
  - c. Pulse: Pulsates for gentle gum stimulation.
- 5 Press the on/off button to turn on the toothbrush.
- 6 Apply light pressure and gently move the brush head slowly across the teeth in a small back and forth motion.
- 7 In order to maximize cleaning effectiveness, divide your mouth in 4 sections. Start with Section 1 and brush each section for 30 seconds. The brushing action will pause between each 30 second interval, indicating when to move on to the next section. The toothbrush will automatically stop after 2 minutes of use. (Fig 3)
- 8 Turn the brush head in direction

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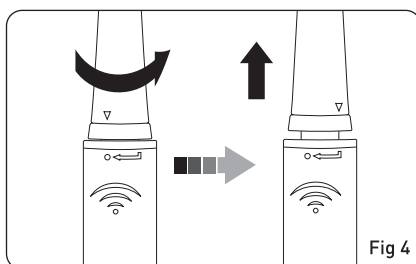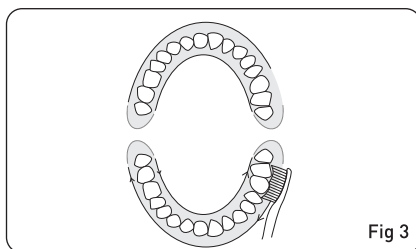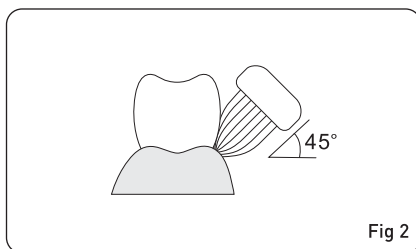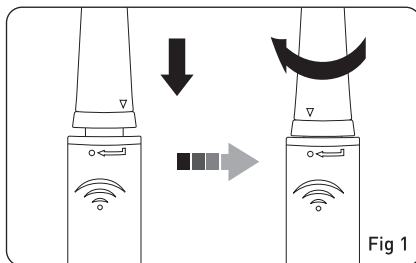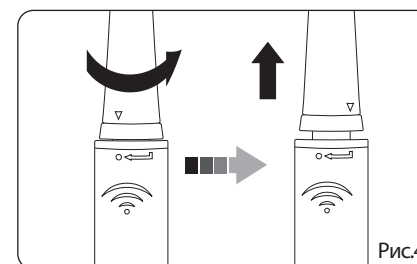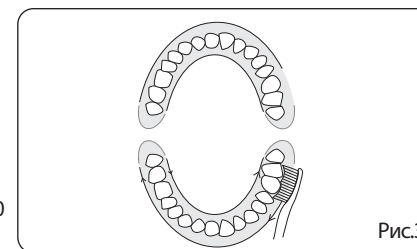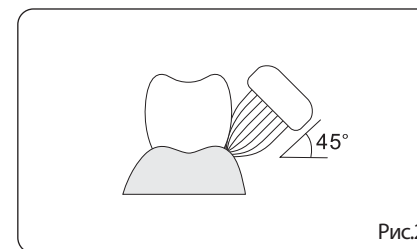
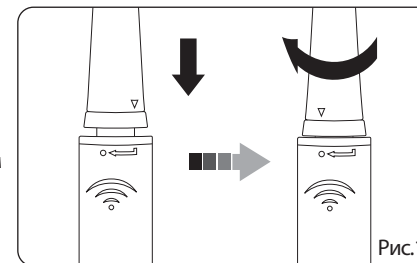
as shown, and then pull it to detach from the toothbrush.  
(Fig 4)


- 9 Rinse both the brush head and handle under running water wipe them dry. (Fig 5)
- 10 Put the handle back into the charger.
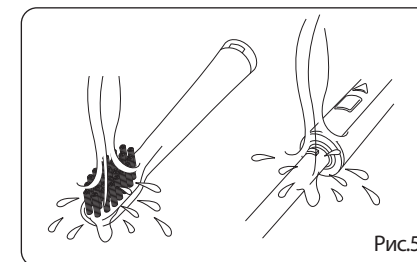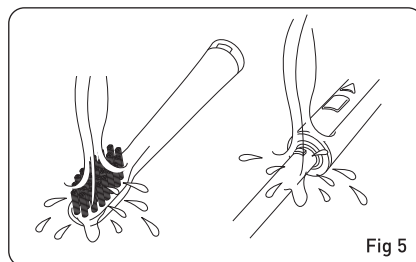
### SANITIZING THE BRUSH HEAD

This product includes a UV sanitizer to clean the surface of your brush heads between use.


- 1 Gently shake the toothbrush to remove any excess water after cleaning.
- 2 Open the UV sanitizer cover.
- 3 Place the brush heads on the stand on the charger top with the bristles directly facing the UV lamp.
- 4 Make sure the charger base with UV sanitizer is plugged into a working outlet.
- 5 Close the UV sanitizer cover to activate the UV sanitizing process.
  - a. The UV sanitizer cover will glow with a slight blue hue indicating the UV light lamp is lit and the sanitizing cycle has

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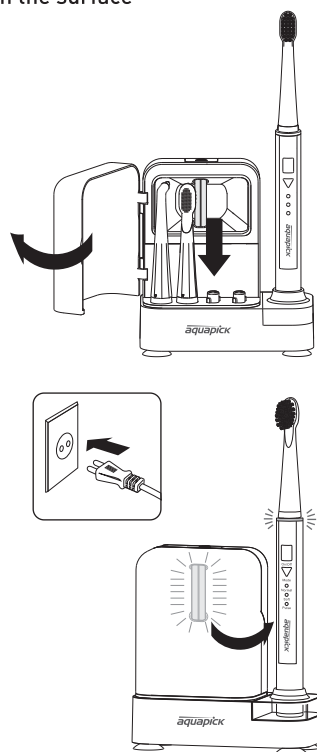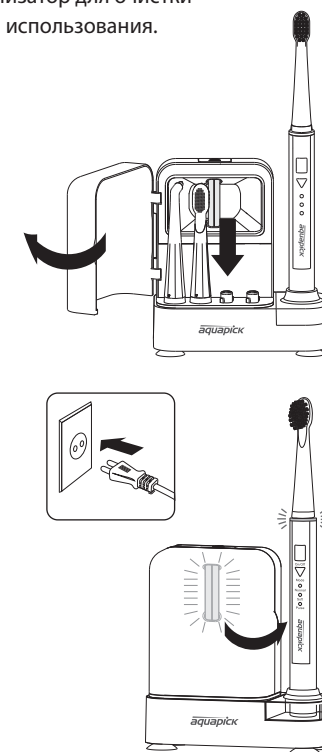- started. Please note this glow may be difficult to view in well-lit areas.
- b. When the sanitizing cycle is complete, the sanitizer will automatically turn off and the cover will no longer glow.
- ⑥ When using the UV sanitizer:
- a. Place only AQ-110 brush heads in the sanitizer.
  - b. The UV sanitizer will run for approximately seven minutes and shut off automatically.
  - c. The UV sanitizer will run automatically every time the UV sanitizer cover is closed.
  - d. Do not place other objects, including brush heads other than AQ-110 into the UV sanitizer.

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
## MAINTENANCE

Be sure to remove the plug from the electrical outlet before cleaning the product.

### Charger

- ① Pull the charger cover to remove.
- ② Using a soft cloth, wipe off any dirt on the charger cover and charger base and UV lamp exterior surfaces.
- ③ Clean the plug periodically to ensure it remains in good working order. The accumulation of dust or dirt on the plug can cause an electric shock, overheater fire.

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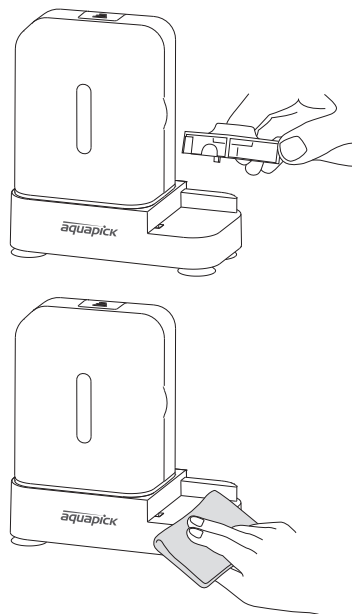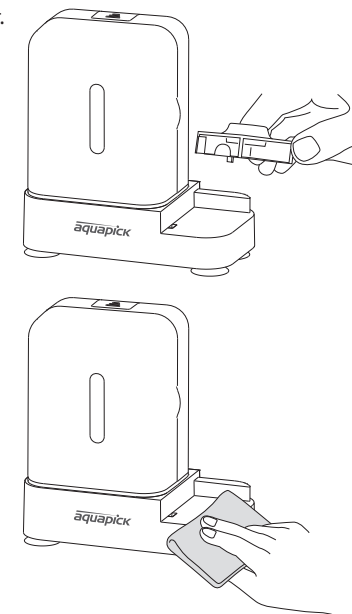
## Toothbrush

- ❶ Remove the brush head.
- ❷ Wipe all surface areas with a damp cloth and place the handle back into the charger. Do not submerge the handle into water.
- ❸ Rinse the brush head and bristles after each use.
- ❹ Do not clean in the dishwasher.

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### TROUBLE SHOOTING

PROBLEM	POSSIBLE CAUSE	SOLUTION
The toothbrush doesn't work.	You have just purchased the product or it has not been used for more than 3 months.	Charge the main unit for at least 24 hours.
Weak or diminished vibration.	The toothbrush is pressed hard against your teeth.	Place the brush head gently against your teeth.
The toothbrush operates only for a short period of time even immediately after being charged.	The battery has reached the end of its life cycle.	Contact service center for support.
	Charging time was insufficient.	Charge the main unit for at least 24 hours.
UV sanitizer light does not turn on.	The sanitizer cover may be open.	Close the sanitizer cover.
	The power cord may not be plugged in.	Plug in the power cord.
	The UV light lamp may be defective.	You may need to replace the UV light lamp.

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## DISPOSAL

The power source of this product is nickel-metal hydride battery which can be recycled. At the end of the product's life, always remove the battery before you discard the product, and recycle or dispose of the battery according to local waste management requirements. Do not throw away the appliance or the UV light lamp (lamp contains mercury) with normal household waste at the end its life. Instead, hand it in at an official waste collection point for recycling. By doing above all, you help to preserve and protect the environment. If you have recycling questions, contact your local waste management office.

### Removing and recycling the battery

Do not disassemble the product at any time other than when disposing of it. Doing so earlier will destroy the product. Observe basic safety precautions. The appliance must be disconnected from the supply mains when removing the battery.

- 1 Before removing the battery, be sure to completely drain the unit of power.
- 2 Use a screwdriver to remove the screw cover and battery door.
- 3 Hold the handle upside down, press the toothbrush shaft, and remove the two hooks to release the internal components of the handle.
- 4 Insert the screwdriver next to the battery connections, and twist to break the connections. Remove the circuit board and pry the battery from the plastic carrier.
- 5 Recycle or dispose of the battery properly. Ask your local waste management office if they accept the remaining toothbrush components for recycling.

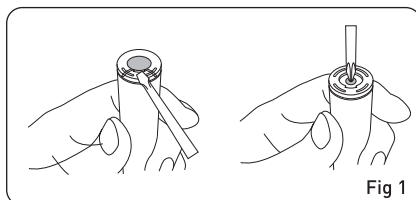
Fig 1

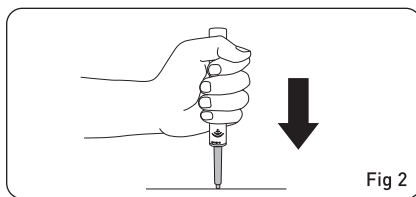
Fig 2

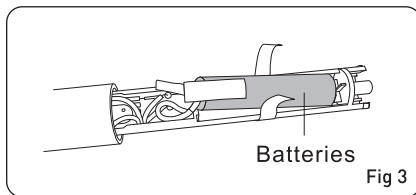
Fig 3

Be careful not to short-circuit the negative and positive poles of the battery. Do not connect the appliance to the mains again after you have removed the battery.

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